

The background features abstract, overlapping green geometric shapes, primarily triangles and polygons, in various shades of green, creating a modern and dynamic visual effect.

Neurograph Project

David Jiashu Wu

Facts about Parkinson's Disease in Australia



What this project is about?

The name tells the answer

Neuro

Graph

Detect Parkinson's
Disease
(Neurology)

By studying
drawing patterns

What this project is about?

Neurograph Project:

Detect the Parkinson's Disease by studying people's drawing pattern.

It's actually a hot topic: From the research papers

present in 25 persons. Furthermore, in two persons (of which one person is familiar with autosomal dominant cerebellar ataxia) drawings were not suited for our post-processing analyses. Therefore, we chose to exclude these results and 1,912 persons were left for further analyses.


FINE MOTOR SKILL ASSESSMENT

Fine motor skill was assessed by requiring participants to trace a picture of a spiral template that was printed on a piece of paper attached to an electronic drawing board (WACOM Graphire Wireless Pen Tablet, model CTE-630BT). Participants were instructed to place the pen in the middle of the spiral before the tracing started (Figure 1A). They were not allowed to lean on the drawing board with their hand or arm. Participants were asked to trace the spiral as accurately and as fast as possible using their dominant hand.

Figure 1.

QUANTITATIVE ANALYSIS OF SPIRAL DRAWING

Automatic quantitative analyses were performed using custom-made software written in MatLab (version 8.1; The Mathworks, Natick, MA, USA). This yielded the following outcome measures: movement time (s), defined by the time it took the participant to trace the spiral; length of drawing (cm), defined as the length of the drawn spiral; *average speed*, defined by the ratio of length of drawing and movement time; *speed variability* (cm/s), defined as the SD of the instantaneous velocity; *deviation from template* (cm²), defined as the area between the template and the drawn spiral; and *number of crossings*, defined as the number of times the drawn spiral crossed the template (Figure 1A). A smoothly drawn spiral with a clinical score of 0A would have a length of drawing about

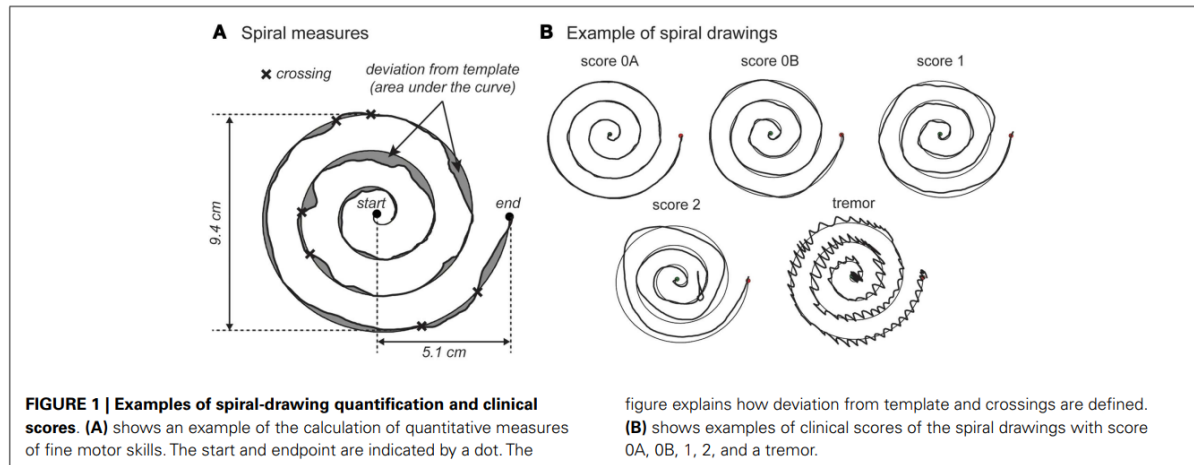

Research paper reference:
<https://www.ncbi.nlm.nih.gov/pmc/articles/PMC4174769/pdf/fnagi-06-00259.pdf>

In the media: Use drawing to detect Parkinson's Disease

BBC Sign in News Sport Weather Shop Earth Travel More Search

NEWS

Home Video World Asia UK Business Tech Science Stories Entertainment & Arts Health World News TV More

Health

Spiral drawing test detects signs of Parkinson's

6 September 2017

f t w e Share



Top Stories

Iran condemns US sanctions move

Tehran dismisses Washington's promise to impose the "strongest sanctions in history".

4 hours ago

Archbishop guilty of child abuse cover-up

1 hour ago

Bongo warns N Korea not to invade

Drawing patterns and Parkinson Detection

What are the connections between them?

Neurograph, Drawing patterns and Parkinson Detection


Parkinson

Shaking, maybe?

Drawing patterns may be potentially influenced

Or just can't focus

The drawing pattern may tell us something interesting

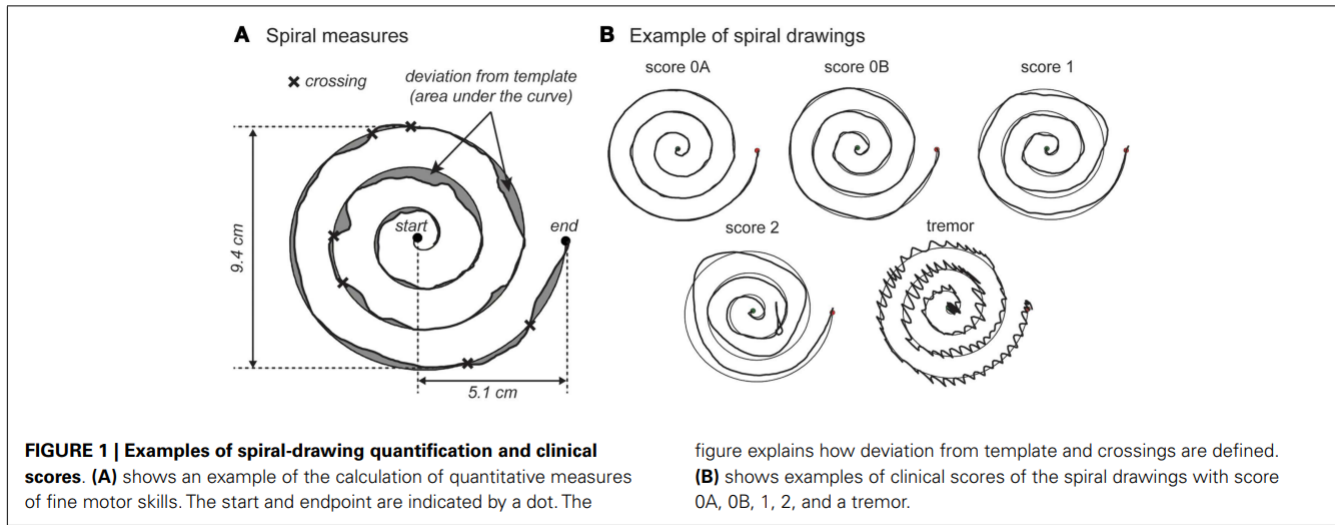

By studying drawing patterns, such as

Average Speed

Drawing Speed changes

Touching pressure

Maximum differences in horizontal direction

Number of crossings

Total drawing time

Deviation from the template image

Analyse them statistically,

We could find some drawing characteristics
that's special for Parkinson patients.

Analyse them statistically,

Then we may be able to distinguish between
normal people and Parkinson patients.

What this project contains

Two parts

Part 1: An Android app which captures the drawing data

Store the data into different format of data files which can be analysed

Part 2: Analyse the drawing data and study interesting features (Done by Cathy)

Introduction to Android app: Neurograph

The background of the slide features abstract, overlapping geometric shapes in various shades of blue, ranging from light sky blue to deep navy blue. These shapes are primarily located on the right side of the slide, creating a modern, tech-oriented aesthetic.

The aim of this app

Capture the data

Store the data, and/or send the data

How the users use it? The Procedure

A new user



Who wants to know whether he is a potential Parkinson patient

Register

Name, age, gender, education level etc should be taken into account

Do the tests

How the users use it? The Procedure

Clinical staff

Output the data by

(1) Sending it via email

(2) Save it locally then connect the device to a PC

How the users use it? The Procedure

Data analyst analyse the data (Cathy)

Then tell whether this user has
Parkinson or not

Done using R and matlab,
Sadly android is not good at data analysis

How the users use it? The Procedure

Maybe after a few months of treatment the user wants to come back to do the tests again


He/she can come back using the registration code
(Something like an user name, it's unique for each user)

Do the tests again to see whether the disease condition changes or not

Drawing tests

In this app we have 6 different kinds of tests, some of them are from the research paper which has been used previously and proved to be useful.


Static Background Test

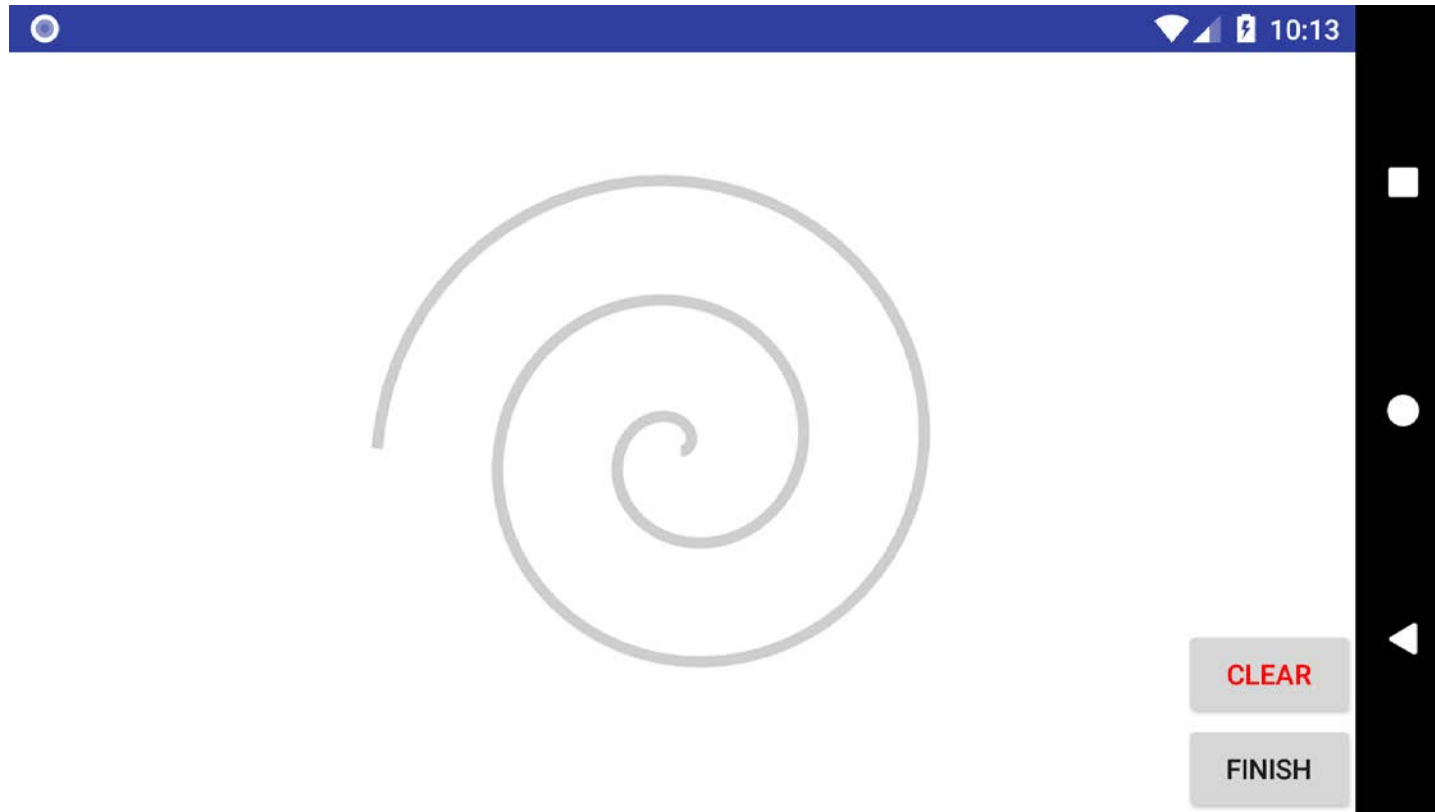


Static Corner Background Test



Dynamic Blank Background Test

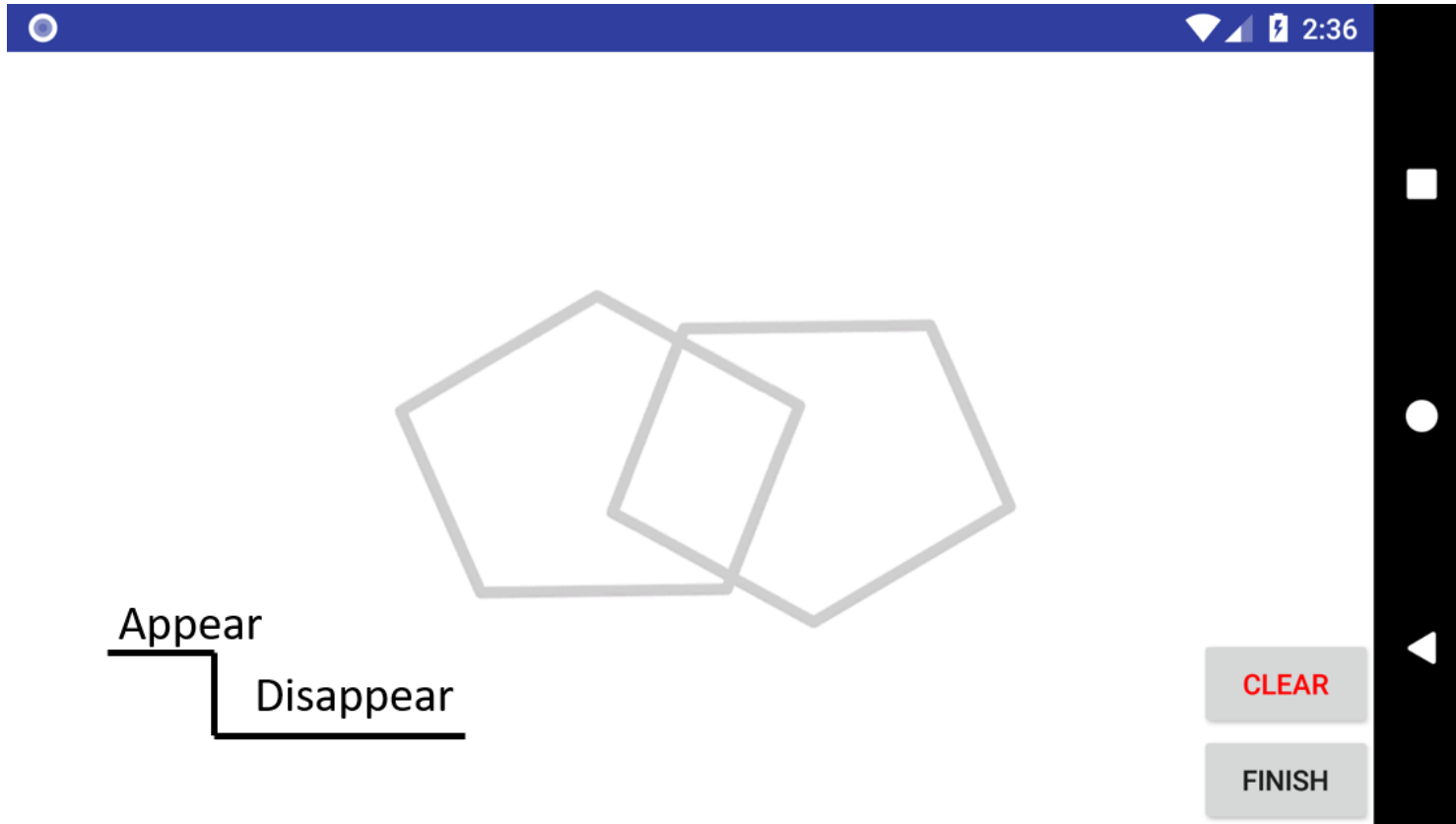


Dynamic Seasonal Background Test

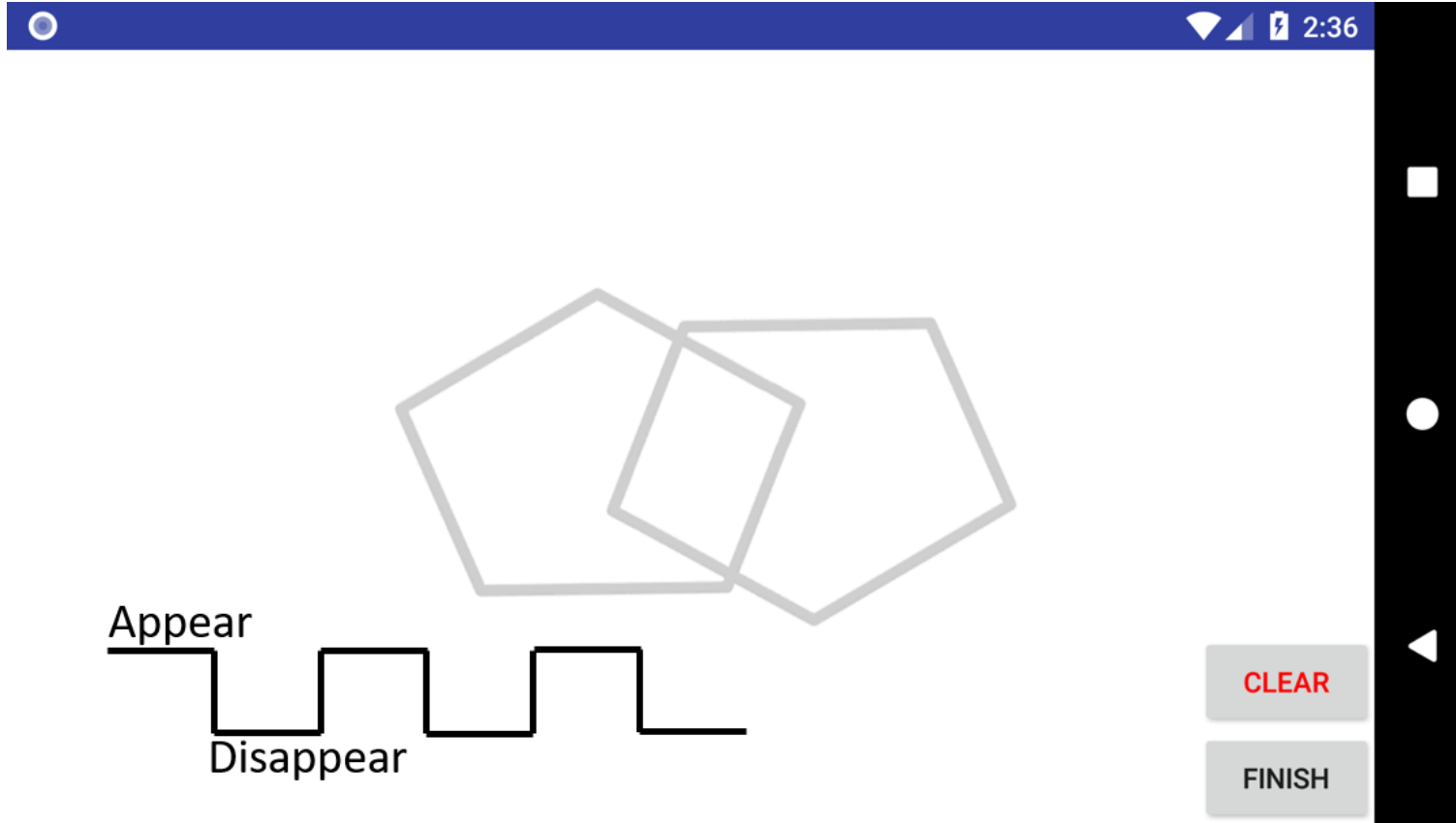

Parallel Line Test



Circular Motion Test

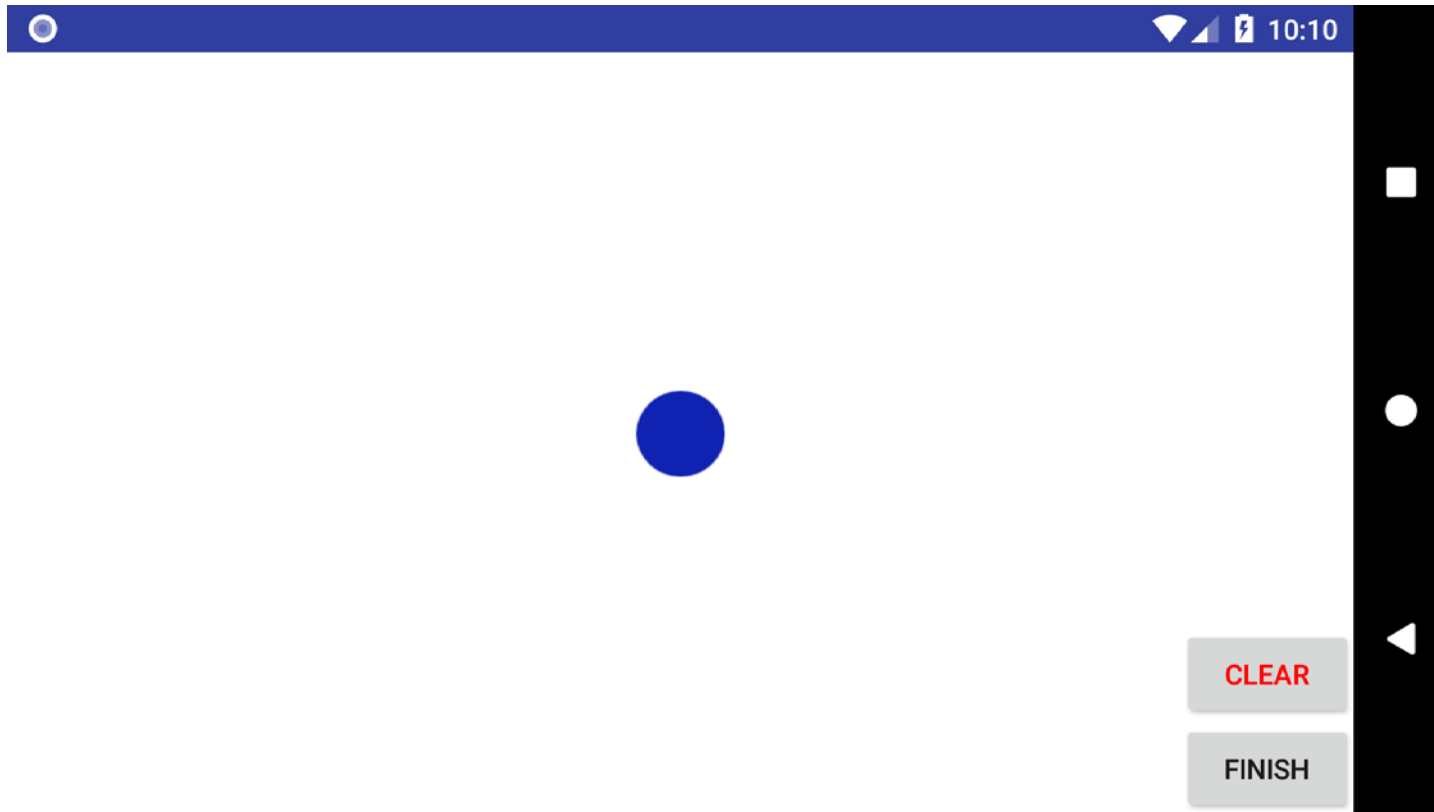

Whether the user can draw the shapes
as required helps us make the decision

Every single detail like deviation area or the changes of speed may tell whether this user has Parkinson's disease or not

Interestingly, android captures dots, not line

In my eye,



In Android's eye,



Actually, Android is smart enough



It captures a point every 0.02 second approximately

The precision can reach 5 decimal places

So there are enough “details” to help us
detect the Parkinson’s Disease

While,


Talk is cheap,
Show me your app...

The app: App demonstration

Collect information about the user

The tests

Output the data

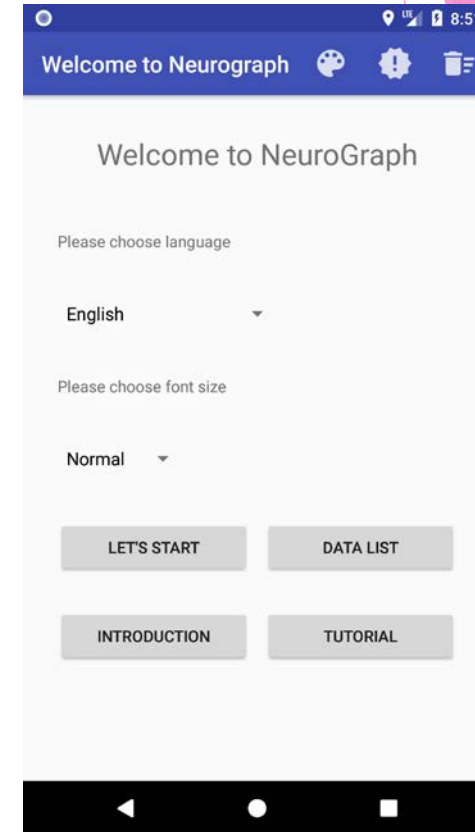

App demonstration

The screenshot shows the Google Play Store interface for the 'Neurograph' app. The top navigation bar includes the Google Play logo, a search bar, and a magnifying glass icon. Below this, the 'Apps' category is selected, with links to 'Categories', 'Home', 'Top Charts', and 'New Releases'. A left sidebar lists 'My apps', 'Shop', 'Games', 'Family', 'Editors' Choice', 'Account', 'My subscriptions', 'Redeem', 'Buy gift card', 'My wishlist', and 'My Play activity'. The main content area displays the app 'Neurograph' by 'MHTI Lab' in the 'Health & Fitness' category. It features a brain icon with 'ng' inside, a 3+ age rating, and a message: 'You don't have any devices.' Below this are 'Add to Wishlist' and 'Install' buttons. At the bottom, four preview images show app screens: 'Welcome to NeuroGraph' (language and font size selection), 'Account Centre' (New User registration and Old User login), 'Full Background Test' (instructions and a 'START TEST 1' button), and 'Pentagon' (a pentagon template).

Google Play

Search

Apps

Categories ▾ Home Top Charts New Releases

My apps

Shop

Games

Family

Editors' Choice

Account

My subscriptions

Redeem

Buy gift card

My wishlist

My Play activity

Neurograph

MHTI Lab Health & Fitness

★★★★★ 2

3+

⚠ You don't have any devices.

Add to Wishlist

Install

Welcome to NeuroGraph

Please choose language

English ▾

Please choose font size

Normal ▾

Account Centre

New User

As a new user, you need to register first.

REGISTER

Old User

Please enter your registration code

Full Background Test

This is the static background (full screen) test. The image will be shown in full screen mode and will appear constantly.

START TEST 1

Corner Background Test

Pentagon

A pentagon template will be

Link to Google PlayStore: <https://play.google.com/store/apps/details?id=com.neurograph.usydjiashuwu.neurograph>

Future TODOs

Identify Signatures

The background features abstract, overlapping geometric shapes in various shades of blue, ranging from light sky blue to deep navy blue. These shapes are primarily located on the right side of the image, creating a modern, layered effect. The rest of the background is a solid, very light blue.

Q & A

Thank you ^_^